

PP & STAINLESS STEEL DRYING RACK



DURABLE • FLEXIBLE • EASY INSTALLATION

PP DRYING RACK

PP Pegboards are injection moulded from polypropylene for structural rigidity, lighter weight and ease of assembly. PP pegboards or drying racks are easy to clean, environmentally friendly and eliminates secondary contamination to hanging glassware.

TECHNICAL

PP Pegboards feature high mobility and design flexibility. All pegs are front side mounted and positioned at a 40 degree upright angle and are removable by lifting upward on the face of the board for cleaning or repositioning of pegs to accommodate larger glassware.

Each pegboard has a pre-moulded drip tray for directing liquids easily and quickly to the drain hose and sink for disposal.

For laboratory space saving and better application, there are two optional pegboards:

- The first is a back to back system for bench mounting
- The second and more universally used is a single system for wall mounting

Each board is accompanied by mounting brackets or connecting brackets for stability for hanging purposes.



PE423-320-52

PP Drying Rack (B) (With 52 Pegs)

Size : 550mm (L) x 125mm (W) x 700mm (H)

Colour : Grey



PE423-310-27

PP Drying Rack (S) (With 27 Pegs)

Size : 400mm (L) x 110mm (W) x 550mm (H)

Colour : Grey

STAINLESS STEEL DRYING RACK

Stainless Steel Pegboards are suitable for use in any laboratory environment. The laboratory pegboards have a durable stainless steel panel and polypropylene pegs with large bases. Our durable construction, allows the laboratory pegboards to last a lifetime while preventing damage to your most delicate glassware.

The stainless steel board can open for cleaning.




PE302-033-27

Stainless Steel Drying Rack With 27 Grey Square Hook Pegs

Size : 420mm (L) x 55mm (W) x 610mm (H)

Colour : Steel

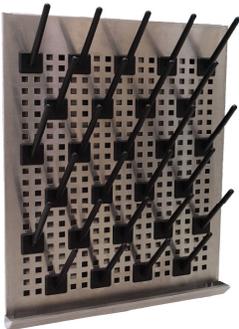

PE302-030-27

Stainless Steel Drying Rack With 27 Black Square Hook Pegs

Size : 420mm (L) x 55mm (W) x 610mm (H)

Colour : Steel



PE302-030

Stainless Steel Drying Rack

Size : 420mm (L) x 55mm (W) x 610mm (H)

Colour : Steel

PP PEG STICK (ROUNDED HOOK)

PP Peg stick with rounded hook is to fit on plastic, phenolic board or wooden board

There are two colour:

- Black
- Grey

**PE121-113**

PP Peg Stick With Rounded Hook (Black)

**PE121-313**

PP Peg Stick With Rounded Hook (Grey)

PP PEG STICK (SQUARED HOOK)



PP Peg stick with squared hook is to fit on steel or stainless steel board


There are two colour:

- Black
- Grey



PE122-113

PP Peg Stick With Squared Hook (Black)

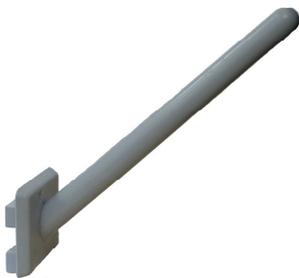

PE122-313

PP Peg Stick With Squared Hook (Grey)